

Multimodal Views of the Human Retina

A. Terry Bahill^{a*} PhD, Patrick J. Barry^b OD

^aSystems Engineering
University of Arizona
1622 W. Montenegro
Tucson AZ USA 85704-1822

^bEnvision Eye Care Center
Tucson Arizona, USA
envisioneyecarecenter@gmail.com

ABSTRACT

Aims: To show the relationships among subject-reported measures of vision, the view through an ophthalmoscope, the view through a slit lamp, a photograph of a face, optical coherence tomography (OCT) images, fundus photos, visual field diagrams and Optomap images.

Methodology: Over 1000 clinical ocular measures (taken on one subject over a six-year period of time) were collected, analyzed and summarized. These measures were reduced to 50 images and tables: they were then categorized and filtered, and the essence resulted in the figures contained in this paper.

Results: This paper shows that the retina literature is full of contradictory nomenclature. For example, clinicians use the term *fovea* to name the 1° diameter disk at the very center of the retina and they use the term *macula* to name the 5° diameter ring that surrounds it. Whereas, anatomists use the term *fovea* to name the 5° diameter disk at the center of the retina and they use the term *macula* to name the 20° diameter ring that surrounds it. This paper demonstrates how the same information appears in the subject's reports of vision, a facial photograph, an optical coherence tomography image, a fundus photo and a visual field diagram. Finally, it shows how to map information between these views.

Conclusions: The retina-viewing techniques analyzed in this paper can be compared qualitatively, but differences in the techniques preclude precise superposition of the images. A perfect mapping is impossible: because (among other reasons) the algorithms for transforming three-dimensional (3D) shapes into two-dimensional (2D) images are nonlinear and are different for different techniques.

Keywords: retina imaging; OCT; fundus photograph,

1. INTRODUCTION

Some ophthalmological diagrams are printed from the subject's point of view and some are printed from the clinician's point of view. For example, the subject would see his left elbow in the lower-left portion of his visual field of view whereas a clinician would see the subject's left elbow in the lower-right corner of a photograph of the subject. Some medical professionals (and their literature) use the term *macula* to name the 5° diameter disk in the center of the retinal, whereas others use the term *fovea* for the same region. Furthermore, there are scores of instruments for viewing the retina. There is no wonder that patients get confused. We think that ophthalmologists and optometrists can use this paper to help explain to their

* Tel.: (520 742-5469 E-mail address: terry@sie.arizona.edu.

27 patients the confusing nomenclature and the overwhelming amount of information generated
28 by their plethora of instruments.

29 Our subject, A. Terry Bahill (ATB), has had a 12 eye surgeries over the last six years. He
30 has been examined and treated by an optometrist, ten ophthalmologists, et al. Over 1000
31 clinical ocular measurements have been performed on his eyes using various retina-
32 visualizing techniques [1]. These measures were analyzed, summarized and reduced to 50
33 images and tables: these data were categorized and filtered: then the essence was put into
34 the figures of this paper. Our objective is to show (in layman's terms) the relationships
35 between the following views: through the subject's eye, through an ophthalmoscope, through
36 a slit lamp, through the clinician's eye looking at a face, fundus photos, optical coherence
37 tomography images, visual field diagrams and Optomap images.

38 1.1 Retinal Distances

39 One millimeter on the human adult retina corresponds to 3.5° of visual angle (0.286 mm/deg)
40 (with a range of 3.4° to 3.6°). The optic disk (the image of the optic nerve and blood vessels
41 entering and leaving the retina) can be modeled as a disk with a diameter of 6.3° (1.8 mm ±
42 0.3). A more complex model is that of an ellipse 6.1° (1.75 mm) wide and 6.5° (1.85 mm)
43 high. On the retina, the center of the optic disk is 15.5° ±1.1° (4.4 mm) nasal and 1.5° ±0.9°
44 (0.4 mm) superior to the center of the fovea [2, 3]. As long ago as 1867 Hulke [4] wrote, "I
45 found the distance of the fovea from the centre of the optic nerve exactly equalled [sic]

46 $1\frac{5}{6}'''$. The triple prime symbol stands for ligne, which equals 2.26 mm. So his measured

47 distance was $1.83''' \times 2.26 \frac{\text{mm}}{'''} \times 3.5 \frac{\text{deg}}{\text{mm}} = 14.5$ degrees. This is within one standard
48 deviation of the distance given above.

49 1.2 Names for Retinal Areas

50 Different classes of people (e.g. clinicians [like optometrists and ophthalmologists],
51 anatomists, engineers, etc.) have named the areas of the retina differently [5]. This has
52 created confusion. Therefore, we created Table 1 and figure 1 to unravel this muddle.

53 **Table 1. Names for areas of the retina, from the center outward to the edge**


Description	Approximate disk diameter	Clinical or Classic name(s)	Anatomical name(s)	Characteristics
Center of foveal pit	0 mm, 0°	Center of fovea	Umbo	
Floor of foveal pit	0.29 mm, 1°	Fovea , fovea centralis [4]	Foveola	It has red and green cones, but no blue cones, no rods and no vasculature
Capillary free zone	0.6 mm, 2°	Foveal avascular zone	Foveal avascular zone	It has no arterioles or venules. It is not circular.
Sloping sides of the foveal pit	1.4 mm, 5°	Macula , macula lutea, yellow spot [4]	Fovea , fovea centralis	It is the darkly pigmented <i>disk</i> in the center of the retina. It is about the size of the optic disk. It has no rods, only cones.
Foveal pit rim-to-rim	1.85 mm, 6.5°	Foveal pit	Foveal pit	The foveal pit rim is the thickest part of the retina.
<i>Ring</i> around	2.9 mm, 10°		Parafovea	Its ganglion cell layer is

the anatomic fovea				composed of four to seven rows of cells.
<i>Ring</i> around parafovea	5.7 mm, 20°		Perifovea	Its ganglion cell layer contains two or three rows of cells.
<i>Disk</i> containing all of the above	5.7 mm, 20°	Posterior pole	Macula , area centralis	
Equator	24 mm, conversion to degrees is meaningless			It is the boundary between the anterior and posterior halves of the eyeball.
The edge of the retina.	30 mm from the center of the fovea.	Ora serrata	Ora serrata	It is on the anterior half of the eyeball. There are no rods or cones anterior to the ora.
For comparison, the optic disk is about 1.8 mm (6.3°) in diameter. For a typical person, the fovea is the size of the fingernail on the little finger with the arm extended.				

54 Column 3 of Table 1 is the classic (or clinical or traditional) retina nomenclature, which goes
55 back at least one hundred and fifty years [4]. Column 4 is the anatomic (or histologic)
56 nomenclature. Dictionaries and encyclopedias usually give the clinical nomenclature.
57 Research papers usually use the anatomic nomenclature. For the last two centuries, the
58 clinical macula (the darkly pigmented 5° disk in the center of the retina) has been called the
59 yellow spot or the macula lutea [4]. With red-free light or certain preparations, it may appear
60 yellow: however, in a fundus photos, it usually looks dark red. Confusingly, a large number of
61 papers describe the macula lutea as a five-degree oval-shaped *yellow spot* in the center of
62 the retina and then they proceed to show a fundus photo with a five-degree dark *red disk*
63 marked as the macula.

64 The foveal pit is a depression in the center of the retina. It was not a regular descriptive
65 feature until cross-sectional views of the retina became common. Now it is a clearly defined
66 landmark. It is easy to measure; the rim-to-rim distance of the foveal pit, 1.85 ± 0.23 mm [6].
67 However, this large standard deviation shows that the inter subject variability of these retinal
68 metrics is large.

69 Figure 1 shows a fundus photo of the retina of the left eye. It has an arc of scar tissue
70 (created by a laser in order to reattach a detached retina) in the upper right corner of the
71 photo and a small area of scar tissue in the lower left corner. The optic disk is circled. The
72 circles in the center show the regions of the retina that are named. The rectangular inset at
73 the bottom is a horizontal cross section of the retina, a slice right through the center of the
74 fovea. On the left side of this rectangle, the surface of the retina starts to dip down into the
75 optic disk.

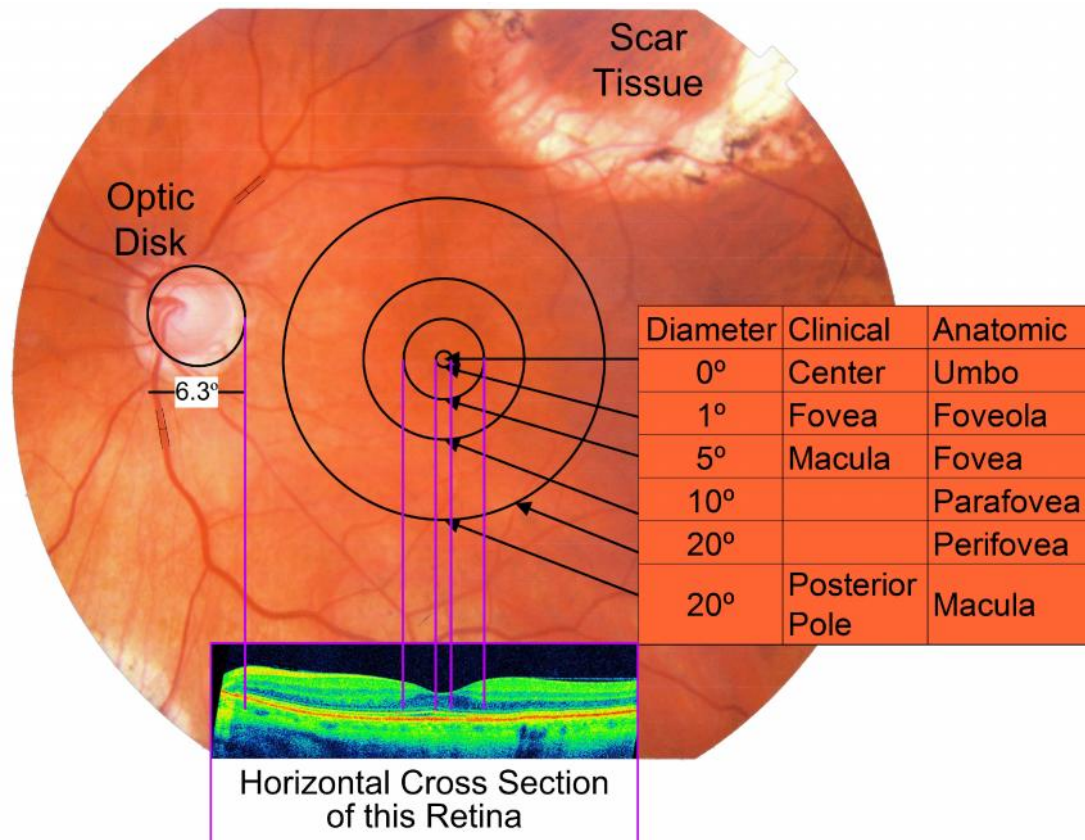

76

77 **Fig. 1. Names for areas of the retina.**

78 **2. METHODS**


79 **2.1 Visualizing the Retina**

80 In this paper, we will consider the following techniques for visualizing the human retina:
 81 subject-reported vision, an ophthalmoscope, a slit lamp, fundus photos, Optical Coherence
 82 Tomography (OCT), visual field diagrams and the Optomap instrument. First, in figure 2, we
 83 show the subject's reported vision of a black-line target composed of concentric squares.
 84 The right eye's image (thin green lines) is the same as the visual target. The left eye's image
 85 (thick red lines) shows the distortion due to a wrinkled retina (macular pucker). The scotoma
 86 in the lower-left quadrant of the figure probably resulted from peeling the inner limiting
 87 membrane (ILM).

88 **2.1.1 Documenting the Wrinkled Retina**

89 Creating figure 2 was complex and time consuming. For this subject, the luminance
 90 threshold for the left eye is about 1.5 log units higher than for the right eye. Therefore, for
 91 these measurements, a 1.5 log unit neutral density filter was fit onto the spectacle lens for
 92 the right eye. The subject would fixate on the center of the target and perceive a small
 93 portion of the target with peripheral vision. Then he would sketch that small portion of figure
 94 2. He would do this with alternate eyes. This process was repeated about one hundred
 95 times. Then a month or so later, the process was repeated again. Details in figure 2 changed
 96 slightly and slowly over the years.

97 The second variant of this process allowed *simultaneous* viewing with the two eyes. Our
 98 subject has strabismus because of the first detached retina operation. So at rest, the image
 99 seen by his left eye drifts to the right. In this measurement, the target was displayed 2.9 m
 100 away from the eyes, so that 5 cm on the target equals one degree on the retina. The neutral
 101 density filter was still on the right lens. Then the subject allowed his eyes to cross so that his
 102 right eye saw the target in its normal position and simultaneously his left eye saw the target
 103 that had drifted to the right. He then the adjusted the sketch of figure 2. The two processes
 104 produced the same sketch for figure 2. This figure is the result of many of measurements
 105 taken over several years. But it was worth the effort, because none of our other retina
 106 visualization techniques gave resolution this fine.

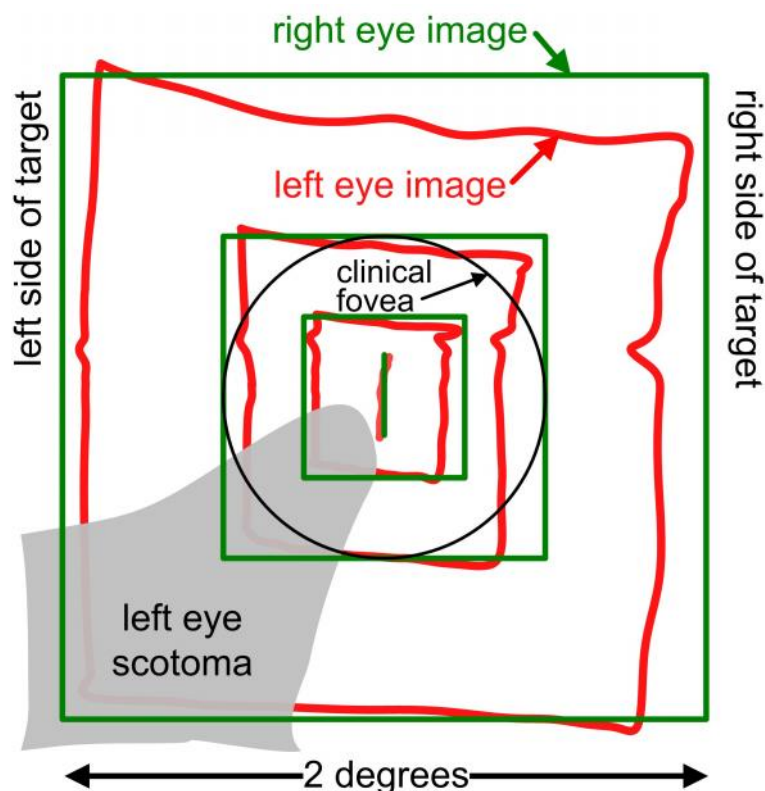

107
 108 **Fig. 2. The subject's viewpoint of a black-line target composed of concentric squares,**
 109 **from [1], reprinted with permission of the authors.**

110 **2.1.2 Techniques for Visualizing the Retina**

111 Next, in figure 3, we compare several other techniques for evaluating the anatomy and
 112 physiology of the retina. The clinician's view of the subject's face is shown in the top row of
 113 figure 3. The image of the subject's right eye (OD) is on the left side of the figure and the
 114 image of the subject's left eye (OS) is on the right side of the figure. With the exception of
 115 this top row in figure 3, all of the data in this paper are for one subject, ATB. Figure 3 shows,
 116 from top to bottom, the clinician's view of the subject's face, OCT images of the right and left
 117 retinas, horizontal cross sections of the retinas, fundus photos of the retinas and visual fields
 118 for the left and right eyes. For each row, the nasal sides of the eyes are in the middle and
 119 the temporal sides of the eyes are on the outside. Each column of figure 3 has the foveae
 120 roughly aligned. The big blue arrows show corresponding locations in the different diagrams.
 121 For the right eye, the blue alignment arrow connects the temporal edges of the optic disk in
 122 the OCT image and in the fundus photo: this is meant to draw attention to the edge of the

123 fovea in the cross section of the retina. For the left eye, this arrow connects, more naturally,
124 the centers of the optic disks. The blue arrows show corresponding locations. Significant
125 points were lined up as well as possible. However, a perfect mapping is impossible: because
126 (among other reasons) the algorithms for transforming three-dimensional (3D) space into
127 two-dimensional (2D) images are nonlinear and are different for different techniques¹.
128 Furthermore, the magnification scales are different for each diagram in this figure.

¹ There are hundreds of algorithms for transforming three-dimensional information (as on the earth or a globe) onto a two dimensional plane (such as a map or an image). One of the oldest and most familiar is the Mercator projection. Google Earth uses a Simple Cylindrical projection. Each of our ophthalmic instruments uses some algorithmic projection, but we do not know which algorithm is used by which medical instrument.

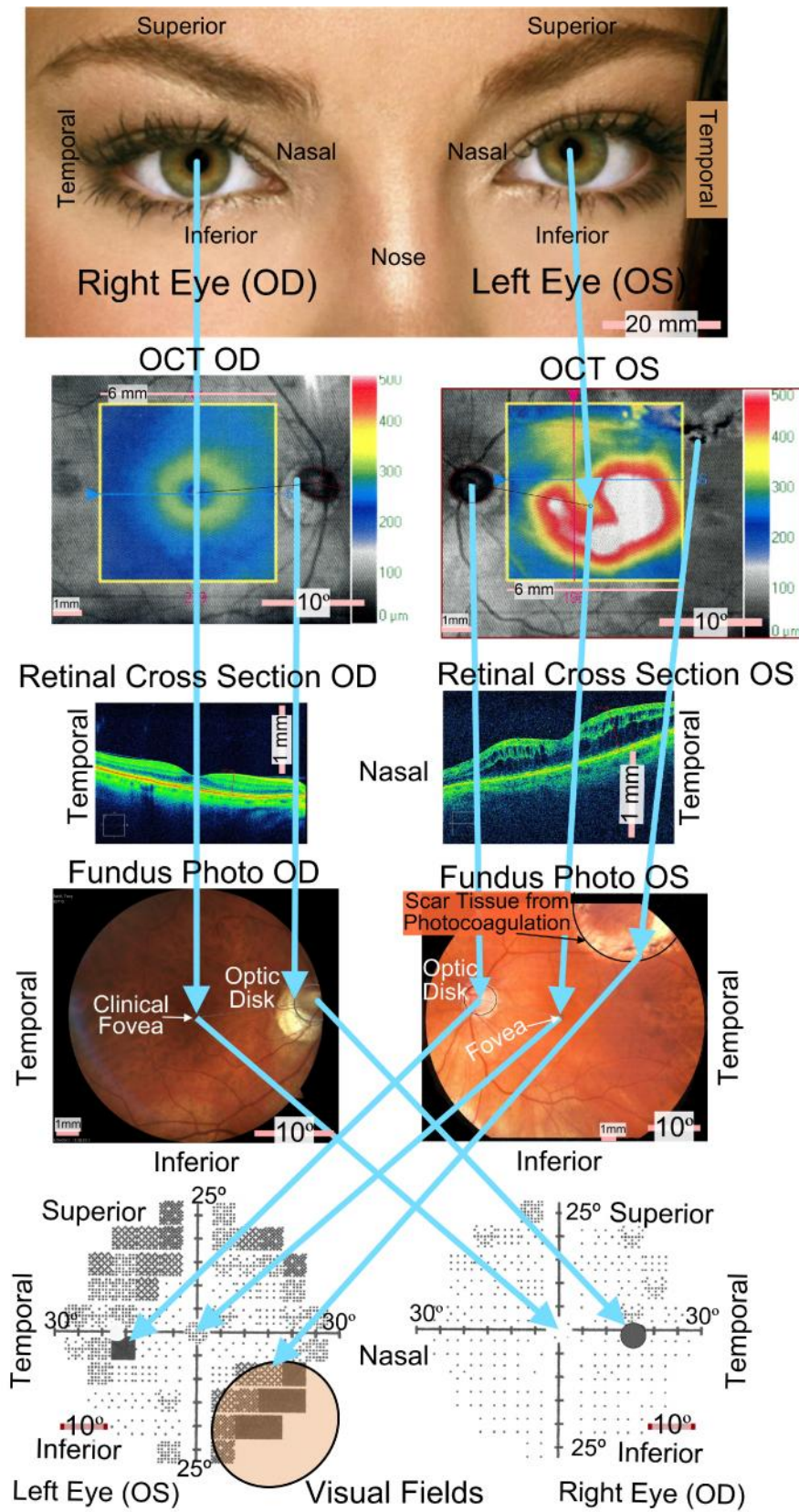


Fig. 3. Shown from top to bottom are the clinician's view of the subject's face, OCT images of the right and left retinas, cross sections of the retinas, fundus photos of the retinas and visual fields for the left and right eyes. Each row has a different scale.

2.2 Fundus Photos

The fourth row of figure 3 shows fundus photos of our subject's right and left retinas. A fundus camera is designed to photograph the interior surface of the eye, including the retina, optic disc, blood vessels, macula and fundus. Using optical lenses, it produces an upright image of the retina for the clinician to see. These images are displayed from the viewpoint of the clinician looking into the subject's face. Therefore, objects that are anatomically in the top part of the retina will appear in the top of the fundus photo. Objects in the nasal part of the retina will appear in the nasal part of the fundus photo. In rows two and four of figure 3, the ugly ring in the superior-temporal region of the left eye's fundus photo and OCT image is scar tissue caused by laser photocoagulation. Because this subject was myopic, the optic disk has a white myopic crescent around it, which is more prominent in this subject's right eye. On fundus photos, arteries are lighter, smaller and less tortuous than veins. The subject does not wear spectacles during fundus photography.

2.3 Visual Field Diagrams

Visual field diagrams show what the subject sees when looking forward out of one eye or the other (monocular viewing). Visual field diagrams are printed from the viewpoint of the subject looking out. This causes a 180° horizontal rotation from the fundus photos. The visual field diagrams of the two eyes are usually plotted side by side, with the left eye on the left of the diagram and the right eye on the right of the diagram, as in the fifth row of figure 3. This puts the temporal parts of the visual fields on the sides of the diagram (near the temples) and the nasal parts in the middle. Things that are in the top part of the visual field will appear in the top of the visual field diagram. Anatomically, the center of the optic nerve (the blind spot) is on average 15.5° nasal and 1.5° superior to the center of the fovea. So visually, the blind spot will be about 15.5° temporal of the center of the fovea as indicated with the solid black circles. The dark squares are areas where the subject does not see. The large dark area in the inferior-nasal region of the left eye's visual field corresponds to the scar tissue in the superior-temporal region of the left eye's fundus photo. In contrast to the other retina-viewing techniques, a visual field plot is linear, because the targets are projected onto a horopter screen (roughly speaking, a section of a sphere.) The subject does not wear his or her glasses when viewing targets in a visual fields instrument; however, an appropriate trial lens should be put in the lens holder in front of the subject's eye.

Visual field diagrams are useful for explaining a subject's vision to the subject. For example, this subject must be on the lookout for things that are above and to his left, like a tree limb above his left shoulder, because now he cannot see the limb with his left eye (because of the scotoma) and he cannot see it with the right eye, because the limb would be out of the field of view of the right eye. He used to be able see objects in this area, but now he cannot. He will not be aware of this deficiency unless it is pointed out to him.

2.4 OCT Images

An Optical Coherence Tomography (OCT) instrument images the retina by measuring the echo time delay and magnitude of reflected light. It directs a beam of light through a beam splitter. described for fundus photos. An OCT instrument images the retina by measuring the echo time delay and magnitude of reflected light. It directs a beam of light through a beam splitter. One beam is focused on the subject's retina, and the other is directed at a reference mirror. Light from the incident beam is reflected off retinal structures at different axial depths and merged with light reflected from the reference mirror. All of this information is put into a three dimensional model in the computer. This model can then be enlarged, rotated,

179 sectioned, etc. The subject does not wear spectacles during OCT measurements: instead
180 the operator focuses the retinal image.

181 The OCT diagrams in the second row of figure 3 show the optic disks (black ovals) and
182 blood vessels radiating from them on the nasal edges. The superimposed central square
183 shows (with a color code) the thickness of the retina (specifically the distance between the
184 inner limiting membrane (ILM) and the retinal pigment epithelium (RPE)). The intersections
185 of the red and blue hairlines in the second row indicate the centers of the foveae, as
186 determined by the computer: for the left eye, this point is wrong, probably because the
187 instrument defines the fovea as the area where the retina is the thinnest and this eye's
188 retinal thickness is altered by edema. Next, row 3 of figure 3 shows a horizontal cross
189 section of the retina passing through the fovea. For the OS there is fluid (macular edema)
190 under the retina: this appears as black columns above the yellow RPE line.

191 We have used three different OCT systems: Heidelberg Engineering, the Optos Optomap
192 and Zeiss Cirrus. The following describes the calibration measures for the Zeiss Cirrus OCT.
193 In a Macular Cube 512 x 128 diagram (second row of figure 3), the ILM-RPE box overlaid on
194 the fundus photo is 6 by 6 mm. The ILM-RPE thickness box (not shown in figure 3) has three
195 concentric circles with diameters of 1, 3 and 6 mm, which correspond respectively to the
196 anatomic fovea, parafoveal ring and perifoveal ring: the anatomic macula disk encompasses
197 all three of these. Finally, on an HD 5 Line Raster (not shown in figure 3), the scan angle,
198 line spacing and length are printed: the default values are zero degrees, 0.25 mm and 6 mm.
199 The correlations between these OCT instruments compare favorably for normal eyes but
200 they differ significantly for abnormal eyes. The two Zeiss instruments, Cirrus and Stratus,
201 had a standard deviation of the mean difference between measurements of central macular
202 subfield thickness of 20 μm [7]. A study of four different models of OCTs devices found an
203 average difference in retinal thickness of 50 μm and a repeatability of 60 μm [8]. Our
204 ophthalmologist feels that a 10% difference between images is noteworthy. This paragraph
205 shows that there is variability within and between OCT instruments made by the same and
206 different manufacturers. Therefore, it takes effort to discover what is being presented, before
207 these results can be compared and contrasted.

208 The point of this section is that similar information about the retina can be obtained from
209 many different instruments. Each instrument has a specific view of the retina. These views
210 are not interchangeable.

211 3. RESULTS

212 3.1 Putting Fundus Photo Information on a Visual Field Diagram

213 A lens flips an image horizontally and inverts it vertically. A fundus photo is taken by a
214 camera looking into the eye and a visual field diagram is derived by the subject looking out
215 of the eye. So to compare a visual field diagram with a fundus photo you must flip the image
216 horizontally and invert it. To flip it horizontally you can move from a position in front of the
217 subject to a position behind the subject. This is the 180° horizontal rotation. However, you
218 did not stand on your head when you went behind the subject, so you must still invert the
219 vertical aspect information.

220 When referring to a region of a photograph or a diagram we will use geometric descriptions,
221 for example the upper-right corner of the photograph. In contrast, when we are referring to
222 anatomy or visual fields, we will use anatomical descriptions, such as the inferior-nasal
223 portion of the left eye's field of view.

224 How would the scar tissue seen in figure 3 affect our subject's visual fields? A large area of
225 scar tissue can be visualized in three images: (1) the upper-right corner of the OS OCT
226 photo, (2) upper-right corner of the OS fundus photo and (3) corresponding lower-right

227 corner of the left eye visual field diagram. Because of the properties of an ocular lens and a
228 fundus camera, scar tissue in the upper-right corner of the fundus photo manifests as a
229 scotoma in the lower-right corner of the visual field diagram. Therefore, when looking
230 straight-ahead (primary position) with his left eye, our subject would not see objects in the
231 inferior-nasal field of view, like a finger tapping his upper lip. Please note that the large white
232 area in the OS OCT image is the area of greatest cystoid macular swelling, it is not the
233 cause of the visual field defect. The cause of the visual field defect is the area of scar tissue
234 that is in the upper-right corner of this image. The other area of scar tissue that is in the
235 inferior-nasal part of the retina is not seen clearly in the photos of figure 3, and therefore it is
236 not discussed here.

237 This paragraph describes the procedure that was used to overlay information from fundus
238 photos and OCT images onto visual field diagrams to produce figure 4. First, we resized one
239 of the images (either the fundus photo or the visual field) so that they had the same scale
240 (degrees of visual angle per millimeter of display). To do this calibration, if we had a clear
241 outline of the optic disk, then we assumed that it was 6.3 degrees in diameter. Otherwise, we
242 assumed that the width of the central retinal artery was $166 \pm 15 \mu\text{m}$ and the width of the
243 central retinal vein was $246 \pm 18 \mu\text{m}$ [9], but these measures had a lot of variability. (Two
244 regions of the central retinal vein are marked in figure 1.) Next, we flipped the fundus photo
245 vertically and put the optic disk from the fundus photo on top of the blind spot in the visual
246 field diagram. Then, to get the rotation correct, we put the fovea from the fundus photo on
247 top of the fixation point (usually the center of the diagram) in the visual field diagram. This
248 process is illustrated in figures 3 and 4.

249 If Photoshop or Visio are being used to illustrate the principle of a lens, then a vertical flip
250 followed by a horizontal flip is equivalent to a 180-degree rotation.

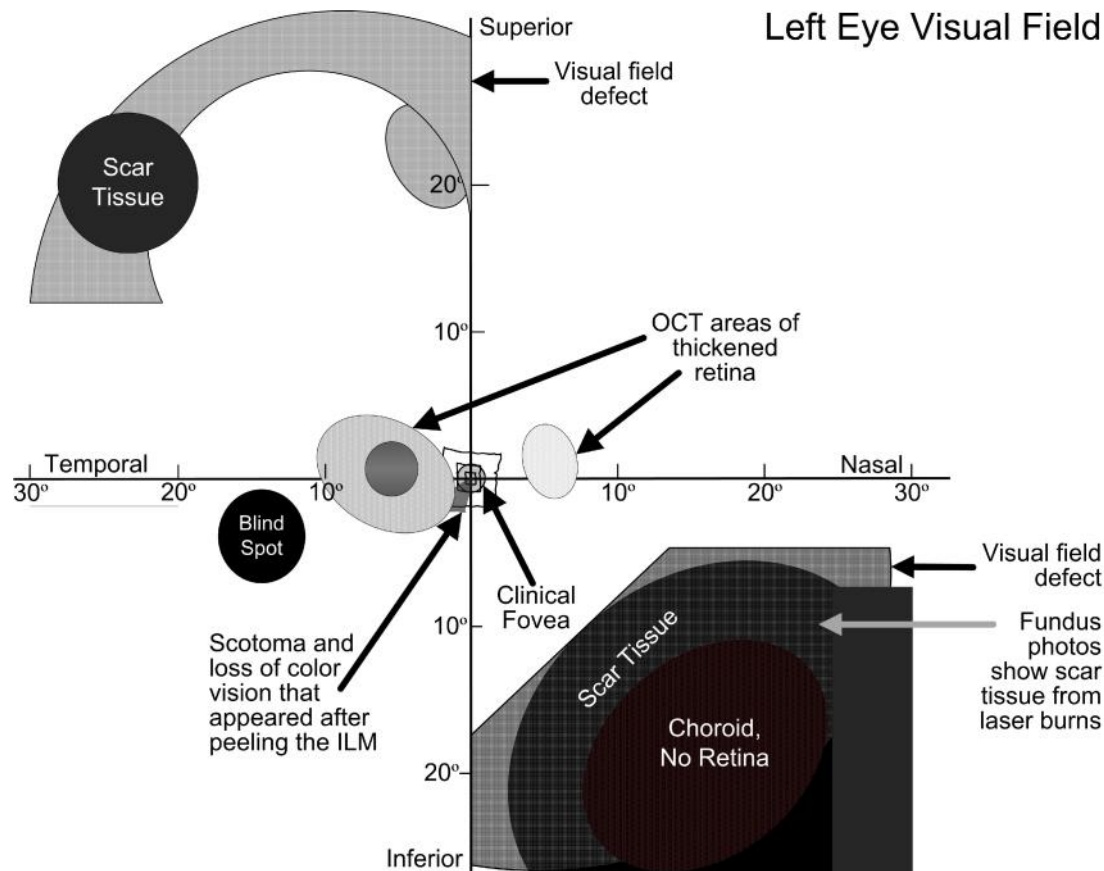




Fig. 4. A visual field diagram for the left eye with fundus photo and OCT information overlaid.

3.2 Optomap Image of the Retina

Figure 5 was designed to show the effects of cyclorotation: it has a traditional fundus photo (top) and an image from a new instrument called the Optomap (bottom), which will be discussed in our Wide-field Imagery section. In both of these images, if we assume that the eye is not cyclorotated, then the optic disk appears to be 15° nasal and 4° superior to the fovea. This vertical distance is larger than the normal $1.5^\circ \pm 0.9^\circ$. Therefore, we re-examined this fundus photo assuming that the eye was cyclorotated by 5°. (Because of his first detached retina operation, our subject has strabismus, which includes 5° excyclorotation [1].) Now, the optic disk appears to be 15.3° nasal and 2.7° superior to the fovea, which is closer to normal. Fundus photos, visual field diagrams and Optomap images all showed this same cyclorotation effect. This effect was not as prominent in OCT images, because the OCT fixation target was a ten-degree asterisk, which gave strong cyclorotation clues and presumably, the subject would suppress his cyclorotational tendency. All of these measurements were made without spectacles, with monocular viewing.

Technical note: Torsion produced by twisting a material in rotational motion is analogous to tension produced by stretching a material in linear motion. The torsional torque is proportional to the angle of rotation just as the tension force is proportional to the distance of

271 extension. So, torsion is proportional to rotation, but torsion is not the same as rotation.
272 However, in the ophthalmology literature cyclotorsion is synonymous with cyclorotation.

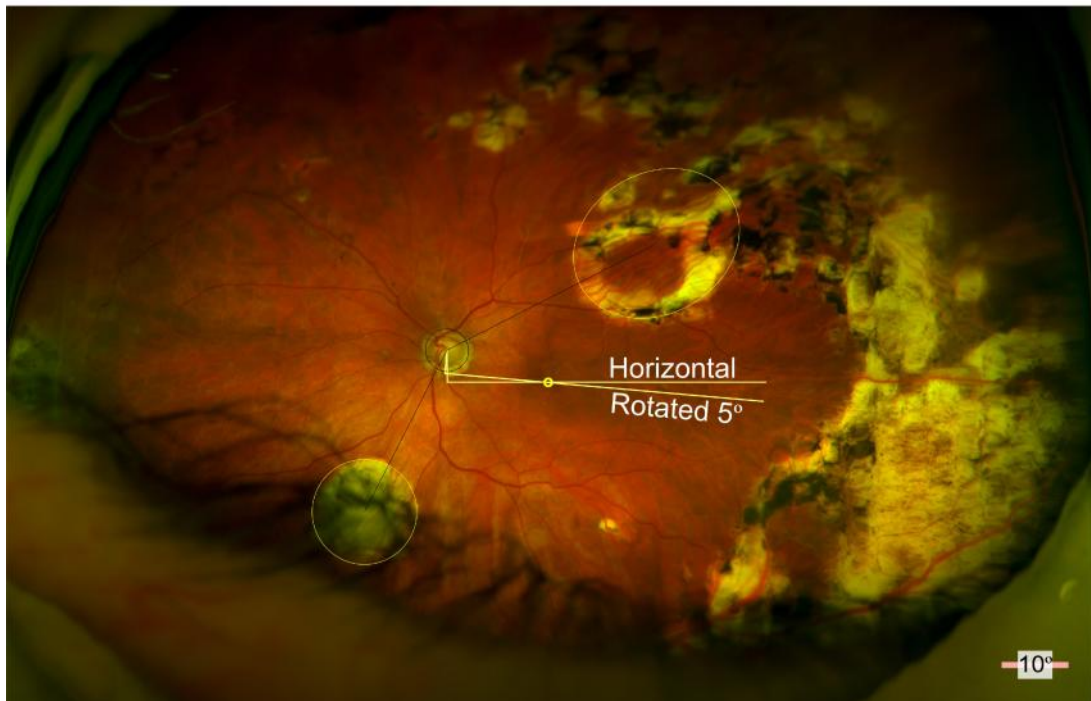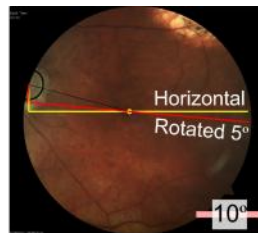

273
274 **Fig. 5. Fundus photo (top) and Optomap image (bottom) of the left eye. Both pictures**
275 **show the effects of cyclorotation of the eye. Circles surround the optic disk, the fovea**
276 **and specific sections of scar tissue.**

277 This section has shown fundus photo information and OCT information being overlaid onto a
278 visual field diagram. If we put this all together, we can show what we believe our subject
279 sees with his left eye as illustrated with figure 4.

280 The macular scotoma of figure 2 cannot be seen in the fundus photos, the OCT images or
281 even the Optomap images. However, in retrospect, we can see hints of this macular
282 scotoma in some visual field diagrams, or perhaps it is just that we see what we believe.

283 3.2 Wide-field Imagery

284 We used a new instrument to image the retina, the Optomap®, which produces images
285 similar to those of a fundus camera. The bottom of figure 5 shows an image computed by
286 this instrument. Compared to fundus photographs, Optomap images (1) have a larger field of
287 view, (2) show features in front of the eye, like eyelashes and eyelids, (3) have greater
288 resolution of arteries and veins and (4) can show arteries and veins in the choroid.
289 Information from these images can be superimposed on visual field diagrams the same as

OCT and fundus photo information. However, the Optomap images are not linear, because they use a nonlinear transformation from a three-dimensional retina to a two-dimensional image. The differences are obvious in figure 5. For both the top and the bottom photos, the size of the calibration bar in the lower-right corner and the size of the optic disk (which is located on the left edge of the fundus photo and in the center of the Optomap photo) are the same. However, the scar tissue on the upper-right margin of the fundus photo is shifted centrally in the Optomap photo. An object that would be “around the backside” of the fundus photo is displayed centrally in the Optomap photo.

The fields of view of our Optomap images were about 25 optic disk diameters horizontally and 13 disk diameters vertically. Assuming the diameter of the optic disk is six degrees, this gives a field of view of 150° by 80°. Three images of normal subjects from the Optos website had fields of view, on average, of 23 disk diameters horizontally and 17 disk diameters vertically, which gives a field of view of 140° by 100°. So we believe that the field of view of the Optomap is around 150° by 90°, which is smaller than the marketing claims.

Figure 5 shows an Optomap image of the left eye and its retina (bottom part of the figure). The unrealistic colors result from using only two scanning lasers: a green laser with a 532 nm wavelength and red laser at 633 nm. The Optomap image shows the optic disk and blood vessels radiating from it. The lower eyelid with its eyelashes is on the top of the image and (because there is an inverting lens between the front of the eye and the retina) the superior part of the retina is also on the top of the image. This left eye image shows scar tissue caused by the laser photocoagulation and cryopexy freezing procedures. In the upper right part of this image, there is a ring of scar tissue about 20° in diameter. Inside of this ring, there is no retina: therefore, the eight wavy red lines are choroid blood vessels. So the Optomap instrument produces photographs that are similar but different from fundus photos. These instruments would be used for different purposes.

3.3 Comparison of Retina-Visualizing Techniques

Table 2 is a cursory comparison of the retina-visualizing techniques that were discussed in this paper. The numbers are subjective. For example, there are many manufactures of fundus cameras: most of them have different fields of view and resolutions. You can use them with or without pupil dilation. Older ones did not produce digital images, some newer ones do. Patents on aspects of fundus cameras run back into the 1920s.

Table 2. Qualitative evaluation of some retina-visualizing techniques

Instrument→ Metric ↓	Subject-reported vision	Ophthalmoscope	Slit lamp	Fundus camera	OCT	Optomap
Invented	a long time ago	1850	1920s	1950s	1990s	2000s
Typical field of view	150° Horz by 135° Vert	20°	2° by 60°	30° to 50°	40°	150° Horz by 80° Vert
Produces digital image	no	no	yes	yes	yes	yes
Produces true color image	yes	yes	yes	yes	no	no
Linear	yes	yes	yes	weakly	claims to be	no
The next four rows use a 0 to 10 scale where 10 is the best.						
Resolution	10	4	6	7	8	9
Usefulness	2	3	6	6	8	9

Patient comfort (without pupil dilation)	10	4	6	7	8	8
Patient comfort (with dilation). Patients dislike dilation and its effects.	not applicable	3	2	4	5	5

The last four rows of Table 2 and all of numbers in Table 3 are subjective. They were established in multiple iterative discussions between the two authors. Barry has all of these instruments in his office and he has made extensive use of each. Other evaluators could easily give different numbers, but we think that the trends would be the same.

Table 3 shows which techniques can evaluate which symptoms. Bahill and Barry [1] give a detailed explanation of each of these symptoms. The right column, doctor's exam, means his or her summary after talking with the patient and reviewing all of the patient's record. The doctor's diagnosis and prognosis could be flawed due to lack of time, lack of documented processes, incorrect test results, failure to notice information, lack of knowledge, poor decision-making, mistakes and different nomenclature. Therefore, the doctor's report could be better or worse than any instrument.

Table 3. Which viewing techniques detect and evaluate which symptoms?

The scale goes from 1 to 10, with a 10 indicating the best detectability for each symptom. Blank cells are zeroes.

Techniques→	Subjects report	Ophthalmoscope	Slit lamp	Fundus camera	OCT	Optomap	Snellen chart	Fluorescein angiography	Visual fields test*	Phoropter	Doctors exam
Symptoms ↓											
Pain in the eye	10										4
Pain in the brain	10										1
Poor visual acuity due to											
wrinkled retina	10 Fig. 2	2	3	4	3	5	2	5	5	7	8
peripheral scar tissue	7	8	5	10 Fig. 3	9 Fig. 3	10 Fig. 5		8	9 Figs. 3 & 4		9
macular scars and holes	9	6	4	8	8	8	7	6	4		7
foveal scotoma	10 Fig. 2						4	2			1
retinal edema	1	3	4	4	10 Fig. 3	5	5	7	5	9	9
macular edema	1	3	3	2	10	6	2	3	3	5	9
color blindness	8										4
Left eye misalignment (phoria)	10									8	9
Double vision (diplopia)	10									8	9
Lack of depth perception (no stereopsis)	9										1

Non-circular and non-responsive (tonic) pupil	5		5							9
Light sensitivity	8	5	5	3						6
*A Macular 10-2 Threshold visual fields test would do better.										

334 In the future, ophthalmological instrument manufacturers will probably add software
335 programs so that their medical instruments can measure more and more of the symptoms in
336 Table 3. For example, the manufacturer could easily add a program to the Optomap so that
337 it could detect non-circular or non-responsive pupils.

338 Table 3 shows that each technique is best for a different symptom: we were surprised at the
339 small overlap. These dozen symptoms seemed like a dozen different problems. But now we
340 can understand the relationships between these problems. We can understand this medical
341 history as one system, instead of a dozen isolated events. Table 4 shows how treating one
342 symptom affects other symptoms. Tables 3 and 4 and figures 3 and 4 brought it all together
343 as one system, instead of a dozen unrelated symptoms. For years, we had a dozen doctors
344 and a dozen symptoms: each doctor treated his most familiar symptoms, one at a time, with
345 his most familiar retina-viewing techniques. Face-to-face communication between the
346 ophthalmologists and optometrists would have led to an earlier understanding of this system.
347 If we had understood this *system* years ago, might the treatment have been different? The
348 patient's eye is a system: it is not a collection of a dozen contradictory symptoms and
349 measures.

350

Table 4. Systems analysis of this case study

The event in this column	prompted the action in this column,	which produced the outcome in this column.	This later ameliorating action	produced this final outcome.
Cataract growth	original cataract operation	pain in the front of the eye	Replacing the intraocular lens (IOL)	eliminated pain in the iris
"	original cataract operation	light sensitivity and a non-circular, non-responsive (tonic) pupil	none is known	permanent damage to iris
Retinal detachment (probably precipitated by the cataract operation.)	many retina operations	macular edema	Peeling the ILM	foveal scotoma, color blindness. wrinkled retina and poor visual acuity.
"	implanting a scleral buckle	inflammation and pain on the back of the eye	Removing exposed sutures on the back of the eye	eliminated pain on the back of the eye
"	implanting a scleral buckle	phoria, diplopia and pain in the brain	Prescribing spectacles with prisms	reduced the pain in the brain
"	implanting a scleral buckle	cyclorotation in rest position causing pain in the brain	none is known	pain in the brain

351 4. CONCLUSION

352 Figure 1 and Table 1 were presented to clarify some very confusing nomenclature that exists
353 in the retina literature: this table has no internal contradictions. Figure 3 and Table 2 showed
354 the relationships between subject-reported vision, a view through an ophthalmoscope, a
355 view through a slit lamp, a photograph of a face, optical coherence tomography (OCT)
356 images, fundus photos, visual field diagrams and Optomap images, all on one unique
357 subject: the subject is unique and the analysis is original. Finally, Tables 3 and 4 showed
358 that the subject's dozen different symptoms were only parts of *one system*.

359 ACKNOWLEDGEMENTS

360 For helpful comments on the manuscript, we thank members of the Systems Engineering
361 Brain Trust: Bruce Gissing, George Dolan, Jerry Swain, Bob Sklar and Brad Sowers.

362 COMPETING INTERESTS

363 Neither author has any conflicts of interest with regards to this paper: furthermore, we
364 believe that there is nothing that could be construed as a conflict of interest. This research
365 was not funded by any organization.

366 AUTHORS' CONTRIBUTIONS

367 Bahill designed the study, performed the statistical analysis, wrote the protocol and wrote the
368 first draft of the manuscript. Barry provided the medical instruments, measurements, clinical
369 images and medical expertise.

370 CONSENT

371 The patient has given his informed consent for this report to be published.

372 REFERENCES

- 373 1. Bahill AT, Barry PJ. A patient's viewpoint of cataract and detached retina surgeries,
374 submitted to *Ophthalmology Review*, 2013;1(2):in press
- 375 2. Rohrschneider K. Determination of the location of the fovea on the fundus, *Investigative*
376 *Ophthalmology & Visual Science*, 2004;45(9):3257-3258.
- 377 3. Kolb H, Nelson R, Fernandez E, Jones BW (editors). *Webvision: The Organization of the*
378 *Retina and Visual System*, Accessed November 2013. Available:
379 <http://webvision.med.utah.edu/book/> .
- 380 4. Hulke JW. On the anatomy of the fovea centralis of the human retina, *Philosophical*
381 *Transactions of the Royal Society of London*, 1867;157:109-115, Accessed November
382 2013. Available: Stable URL: <http://www.jstor.org/stable/108970> .
- 383 5. Park SS. The anatomy and cell biology of the retina, in *Duane's Foundations of Clinical*
384 *Ophthalmology*, Volume 1, Chapter 19, Lippincott Williams & Wilkins, 2006, Accessed
385 November 2013. Available:
386 <http://www.oculist.net/downatn502/prof/ebook/duanes/pages/v7/v7c019.html> .
- 387 6. DuBois AM, Hansen BR, Cooper RF, Beringer J, Dubra A, Carroll J. Relationship
388 between the foveal avascular zone and foveal pit morphology, *Investigative*
389 *Ophthalmology & Visual Science*, 2012;53(3):1628-1636.

- 390 7. Lammer J, Scholda C, Prünte C, Benesch T, Schmidt-Erfurth U, Bolz M. Retinal
391 thickness and volume measurements in diabetic macular edema: a comparison of four
392 optical coherence tomography systems. *Retina*. 2011;31(1):48-55.
- 393 8. Giammaria D, Ioni A, Bartoli B, Cofini V, Pellegrini G, Giannotti B. Comparison of
394 macular thickness measurements between time-domain and spectral-domain optical
395 coherence tomographies in eyes with and without macular abnormalities. *Retina*.
396 2011;31(4):707-716.
- 397 9. Taarnhøj N, Larsen M, Sander B, Kyvik KO, Kessel L, Hougaard JL, Sørensen T.
398 Heritability of retinal vessel diameters and blood pressure: a twin study, *Invest.*
399 *Ophthalmol. Vis. Sci.* 2006;47(8):3539-3544

400 ACRONYMS

401	HD	High Definition
402	ILM	Inner Limiting Membrane
403	OCT	Optical Coherence Tomography
404	OS	Left eye
405	OD	Right eye
406	RPE	Retinal Pigment Epithelium

407 BIOGRAPHICAL SKETCHES

408 **Terry Bahill** is an Emeritus Professor of Systems Engineering and of Biomedical
409 Engineering at the University of Arizona in Tucson. He received his Ph.D. in electrical
410 engineering and computer science from the University of California, Berkeley, in 1975. Bahill
411 has worked with dozens of high-tech companies presenting seminars on Systems
412 Engineering, working on system development teams and helping them to describe their
413 Systems Engineering processes. He holds a U.S. patent for the Bat Chooser, a system that
414 computes the Ideal Bat Weight for individual baseball and softball batters. He was elected to
415 the Omega Alpha Association: the Systems Engineering Honor Society. He received the
416 Sandia National Laboratories Gold President's Quality Award. He is a Fellow of the Institute
417 of Electrical and Electronics Engineers (IEEE), of Raytheon Missile Systems, of the
418 International Council on Systems Engineering (INCOSE) and of the American Association
419 for the Advancement of Science (AAAS). He is the Founding Chair Emeritus of the INCOSE
420 Fellows Selection Committee. His picture is in the Baseball Hall of Fame's exhibition
421 "Baseball as America." You can view this picture at <http://www.sie.arizona.edu/sysengr/>.

422 **Patrick J. Barry** earned his BS at Siena College in Loudenville, NY, with a major in biology
423 and a minor in chemistry. He earned a Doctor of Optometry Degree from the Pennsylvania
424 College of Optometry in 1977 and has been in private practice since 1989, providing primary
425 care Optometric services to the Tucson community. Dr. Barry is a member of the American
426 Optometric Association, the Arizona Optometric Association and the Southern Arizona
427 Optometric Society. He is a Past- President of both the Southern Arizona Optometric Society
428 and the Arizona Optometric Association, where he received the prestigious Optometrist of
429 the Year Award for his dedication to the Optometric profession.
430 envisioneyecarecenter@gmail.com

431